

# clinell®

[www.clinell.com](http://www.clinell.com)

**GOOD  
DESIGN**

## Only 10 pieces

which quickly and easily assemble  
onto the durable, stainless steel frame -  
containing no screws, nuts or bolts

Our award-winning\* commode features smooth, moulded, non-porous plastic parts and unique EasyMove™ and TotalStop™ lockable wheels. The patented 360° rotating and removable footrests incorporate a patented 'Clip & Slide' mechanism. The latex-free seat cover contains no stitching or foam components which are known to harbour microorganisms.

This minimalist design reflects the latest thinking in infection prevention, making it easy to clean, which in turn increases compliance and has been used by over 70 hospitals in the UK in the last year alone.


\* Winner of the Good Design Awards 2012,  
organised by The Chicago Athenaeum  
Museum of Architecture and Design



# EASYS CLEAN COMMODE

T: 020 7993 0030

E: info@clinell.com

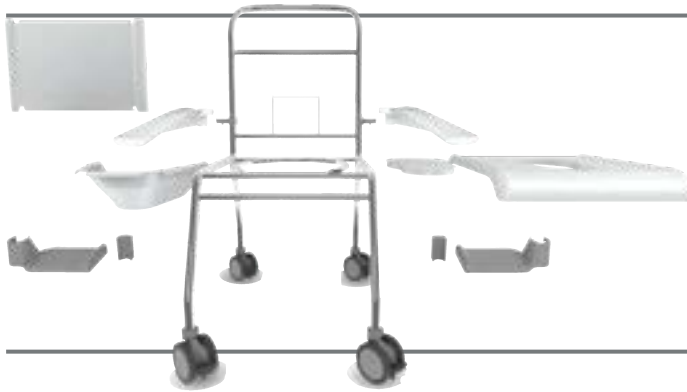

Our unique, super hygienic, multi-functional commode. It can be fully dismantled, making it simple and easy to clean.

## Good Design Award

Due to its unique, superior design features, Clinell's EasyClean Commode has been awarded with the prestigious GOOD DESIGN™ award - the world's oldest, most recognised, annual Design Awards by The Chicago Athenaeum Museum of Architecture and Design.

## Multi-functional

Suitable for both pulp and plastic bedpans. Can also be used as a shower chair or over a toilet.

## Removable single piece, plastic bedpan rack

The fully removable rack allows for easy cleaning, over the toilet access and the insertion of top-loading pulp bedpans.

## Durable, stainless steel frame

3 year frame warranty.

## Easy move lockable wheels

With patented EasyMove™ and TotalStop™ design\*.

## Hinged arm rests

Enables easy patient transferring and support handling. The arms lock securely when down.

## Non-stain plastic

Resistant to staining and chemical damage.

## Easy clean

No screws, nuts or bolts and smooth, curved, non-porous surfaces allow for better and easier cleaning which increases compliance. This means fewer places for dirt and germs to accumulate.

## Rotating, removable footrests

Rotate 360° so they can be turned out of the way for easy mounting or dismounting. Incorporates patented Clip & Slide mechanism for easy removal\*\*.

## Safe and sanitary plastic seat

Latex-free seat cover, which sits flush within the seat, with no stitching or foam which guards against germ build up, allergies and is easy to clean.

## Convenient seat storage

The seat cover can be neatly stored at the back of the commode.



Plastic or pulp bed pans slide in from the rear.



Top loading, moulded, single piece, plastic bedpan rack.



Hinged, lockable arms for easy patient transfer.



Seat cover stowed safely at the back of the seat.



Patented EasyMove™ and TotalStop™ design wheels.



Patented clip and slide footrests.

SPECIFICATION	SIZE
Height	105cm (41.3")
Width	60cm (23.6")
Depth	77cm (30.3")

SPECIFICATION	SIZE
Width between arms	49cm (19.3")
Seat clearance to floor	51cm (20")
Maximum user weight	190kg (30st)

PRODUCT	UNIT OF ISSUE	ORDER CODE
EasyClean Commode	Single unit	CCOM1

\* Function patent: ZL 2006 2 0062320.0 - Appearance patent: ZL 2006 3 0066473.8

\*\* Design patent: GB 2496861B

JBN18714